


CLICK-IN SYSTEM



MULTIFUNCTIONAL & FLEXIBLE

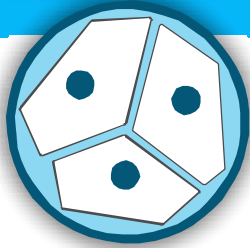

Innovative system for the fixation of cell culture vessels on xy-stages and scanning stages at inverse microscopes

APPLICABLE WITH

Universal Mounting Frame KM
010-800 486

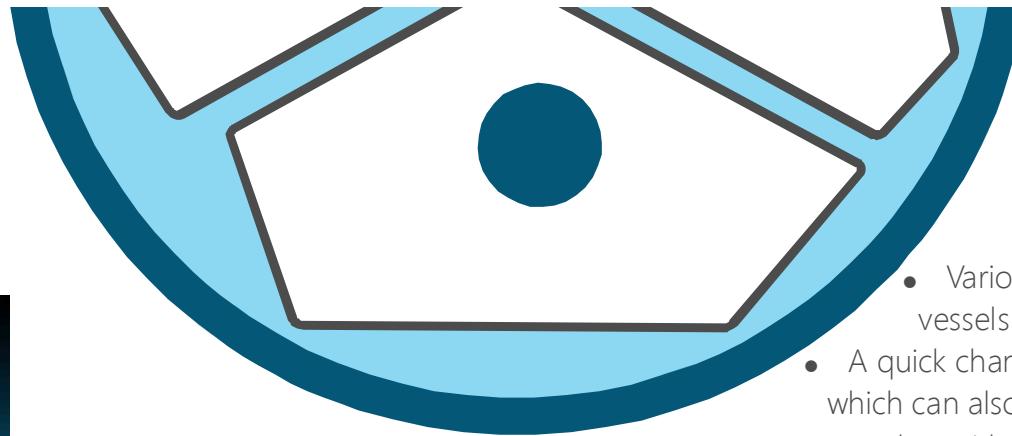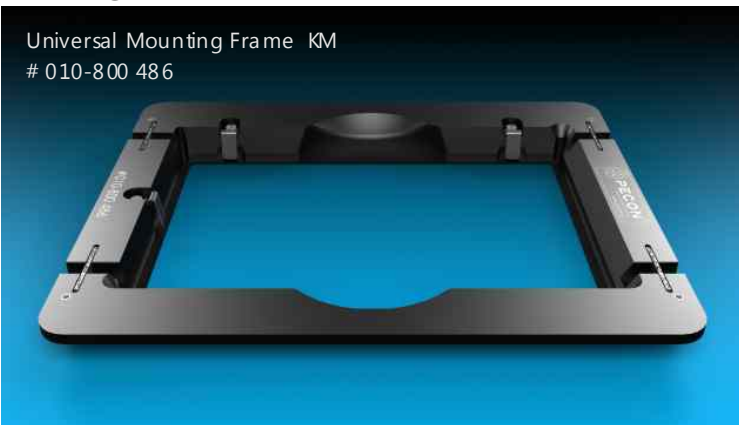

There is a number of advantages in comparison to the conventional Universal Mounting Frames:

- Various different types of cell culture vessels can be used.
- A quick change of the cell culture vessels, which can also be inserted and taken out together with the insert plate, is possible.
- A firm position of the cell culture vessel by clamping springs in the frame and in the Z-direction by means of spring clips (easy to assemble or disassemble) is ensured.
- Cost-efficient alternative to other Universal Mounting Frames.
- With a CO₂-Cover and a CO₂-Controller, the pH-value in the nutrition medium can be controlled.
- Custom-specific designs are possible anytime.

PeCon's innovative **Click-In System** combines a wide range of application with a simple handling. It provides for a firm fixation and a stable position of the cell culture vessels during the observation under the microscope.



CELL CULTURE VESSEL

INSERT PLATE

UNIVERSAL MOUNTING FRAME

READY-TO-USE CLICK-IN SYSTEM

The Universal Mounting Frame KM serves as a base frame (suitable for xy-stages and scanning stages with a cut-out of 160 x 110 mm) into which multiplates and different insert plates can be clicked in. Furthermore, custom-specific designs or modifications are possible anytime.

PECON -
SYSTEMS FOR
CELL CULTURE
AND MICROSCOPY

INSERT PLATES **CLICK-IN SYSTEM**

CLICK-IN SYSTEM PREMIUM

The **Click-In System Premium** consists of a 6 mm thick black eloxadized aluminium plate. The different cell culture vessels are laterally fixed by clamping springs and can be additionally fixed in the Z-direction from above by spring clips.



35er Petri dishes
(Ø 30 - 40 mm)

000 440



"60" Petri dishes
(Ø 47 - 56 mm)

000 415



2x "35" Petri dishes
(Ø 30 - 40 mm)

000 441



1x "35" Petri dish
(Ø 30 - 40 mm) &
1x "60" Petri dish
(Ø 47 - 56 mm)

000 427



POC-R2 &
POCmini-2
by LaCon

000 442



all object slides
and chamber sys-
tems of different
manufacturers

000 470



"100" Petri dishes
(Ø 78 - 90 mm)

000 482

FURTHER INFORMATION:

PeCon GmbH
Ziegeleistrasse 50
89155 Erbach
Germany

Email: info@pecon.biz

Web: www.pecon.biz