
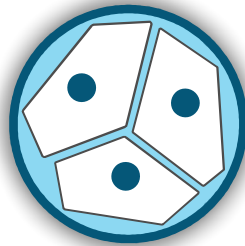


# CLICK-IN SYSTEM



## MULTIFUNCTIONAL & FLEXIBLE

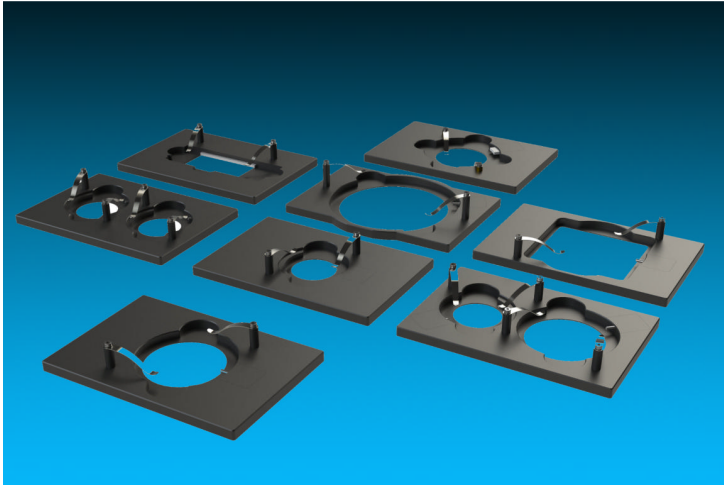

Innovative system for the fixation of cell culture vessels on xy-stages and scanning stages at inverse microscopes

### APPLICABLE WITH

Universal Mounting Frame KM  
# 010-800 486

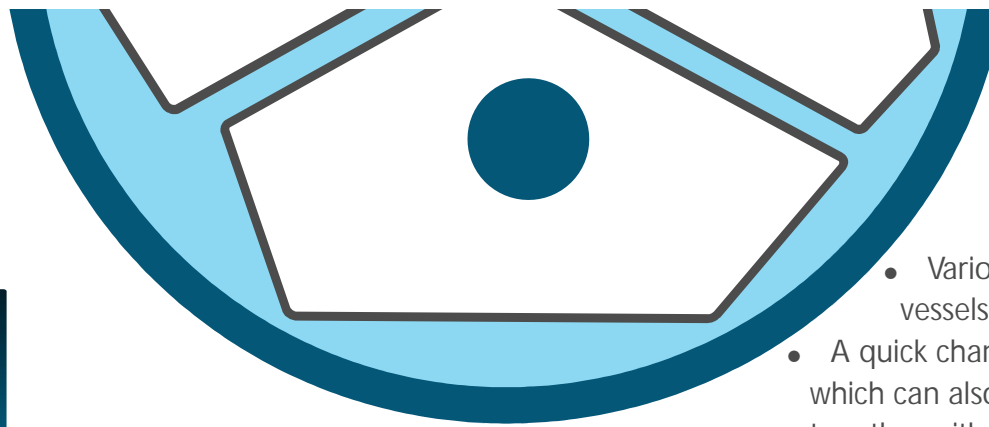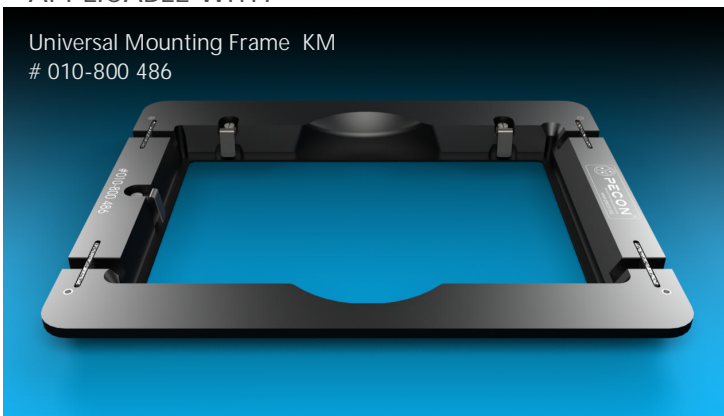

PeCon's innovative Click-In System combines a wide range of application with a simple handling. It provides for a firm fixation and a stable position of the cell culture vessels during the observation under the microscope.



CELL CULTURE VESSEL

INSERT PLATE

UNIVERSAL MOUNTING FRAME

READY-TO-USE CLICK-IN SYSTEM

The Universal Mounting Frame KM serves as a base frame (suitable for xy-stages and scanning stages with a cut-out of 160 x 110 mm) into which multiplates and different insert plates can be clicked in. Furthermore, custom-specific designs or modifications are possible anytime.

There is a number of advantages in comparison to the conventional Universal Mounting Frames:


- Various different types of cell culture vessels can be used.
- A quick change of the cell culture vessels, which can also be inserted and taken out together with the insert plate, is possible.
- A firm position of the cell culture vessel by clamping springs in the frame and in the Z-direction by means of spring clips (easy to assemble or disassemble) is ensured.
- Cost-efficient alternative to other Universal Mounting Frames.
- With a CO<sub>2</sub>-Cover and a CO<sub>2</sub>-Controller, the pH-value in the nutrition medium can be controlled.
- Custom-specific designs are possible anytime.

PECON -  
SYSTEMS FOR  
CELL CULTURE  
AND MICROSCOPY

# INSERT PLATES CLICK-IN SYSTEM


## CLICK-IN SYSTEM PREMIUM

The Click-In System Premium consists of a 6 mm thick black eloxadized aluminium plate. The different cell culture vessels are laterally fixed by clamping springs and can be additionally fixed in the Z-direction from above by spring clips.

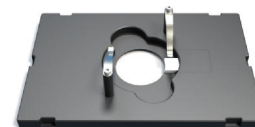

35er Petri dishes  
(Ø 30 - 40 mm)

# 000 440

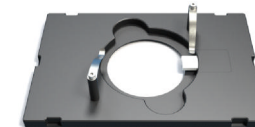

"60" Petri dishes  
(Ø 47 - 56 mm)

# 000 415




2x "35" Petri dishes  
(Ø 30 - 40 mm)

# 000 441



1x "35" Petri dish  
(Ø 30 - 40 mm) &  
1x "60" Petri dish  
(Ø 47 - 56 mm)

# 000 427

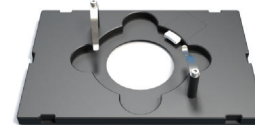

POC-R2 &  
POCmini-2  
by LaCon

# 000 442



Insert Plate  
for 4-Well  
Multiwell plates

# 001 236



"100" Petri dishes  
(Ø 78 - 90 mm)

# 000 482



All object slides  
and chamber sys-  
tems of different  
manufacturers

# 000 470

### FURTHER INFORMATION:

PeCon GmbH  
Ziegeleistrasse 50  
89155 Erbach  
Germany  
Email: [info@pecon.biz](mailto:info@pecon.biz)  
Web: [www.pecon.biz](http://www.pecon.biz)