
PECON®
Systems for cell culture and microscopy

"Open" perfusion of cell cultures

Perfusion-Set for "60" Petri dishes
500-800260

Perfusion-Set for "60" Petri dishes

consisting of:
cover with glass insert, # 000861; perfusion ring, # 501466;
canal tube inlet, # 501425; canal tube outlet, # 501426;
grub screws # 900094; screwdriver # 900405

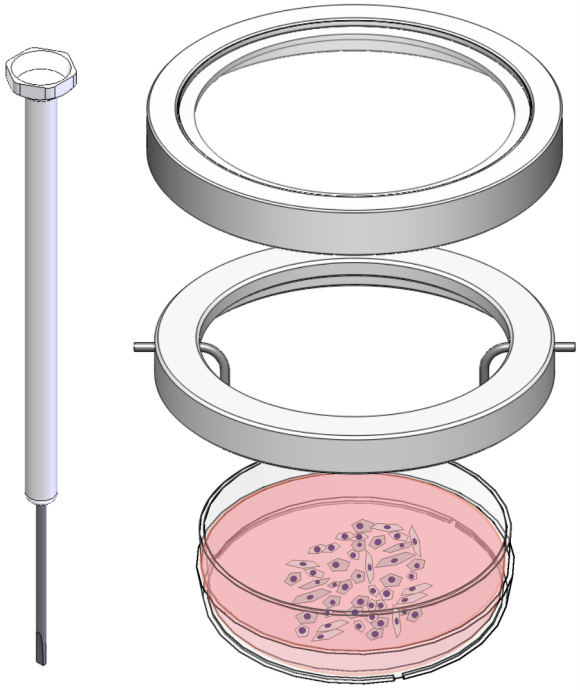cover with glass insert

- stainless steel V4A
- sterilizable
- glass insert
- cut-outs for the canal tubes

perfusion ring


- stainless steel
- sterilizable
- lateral openings for the canal tubes
- Screws for the fixation of the canal tubes

Petri dishes "60"

- different manufacturer with outer diameter 50.5 - 56.0 mm

Scheme of an "open perfusion"

assembly and procedure of the perfusion

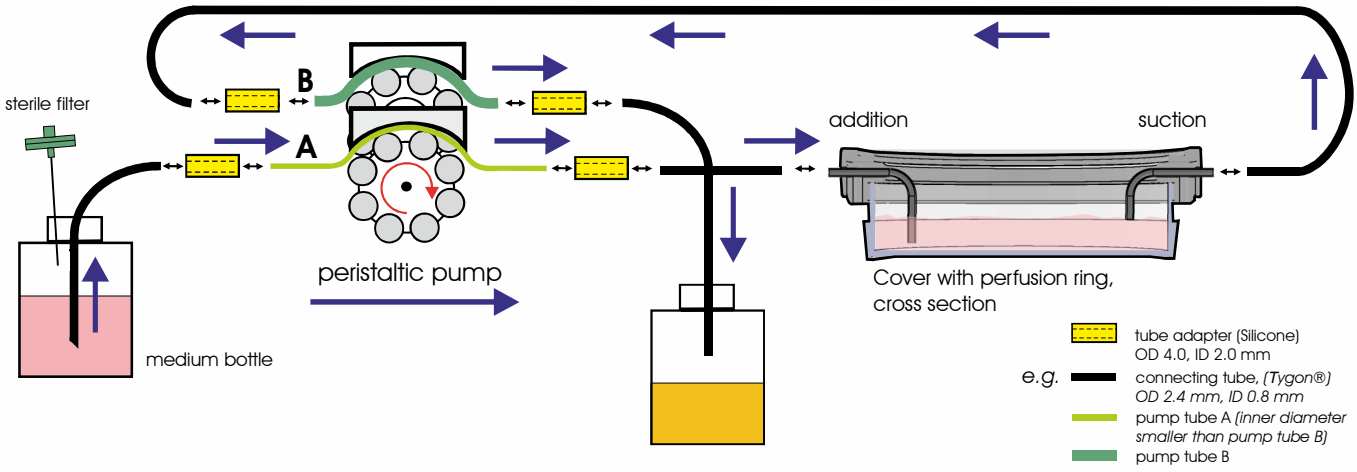This Perfusion-Set is compatible to a large variety "60" Petri dishes.
Please do not hesitate to contact us regarding a cover with perfusion adapter for other Petri dishes!

Version 1.2, 03-2019
www.pecon.biz


PECON®
Systems for cell culture and microscopy

"Open" perfusion of cell cultures


Perfusion-Set for "60" Petri dishes
500-800260

Specifications


cover with glass insert

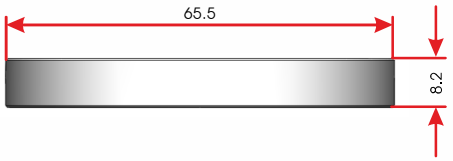perfusion ring with canal tubes

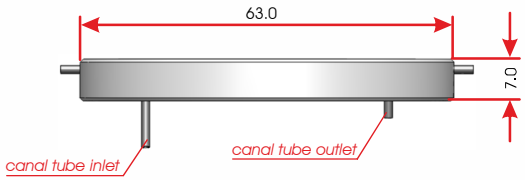
perfusion ring with cover - bottom view

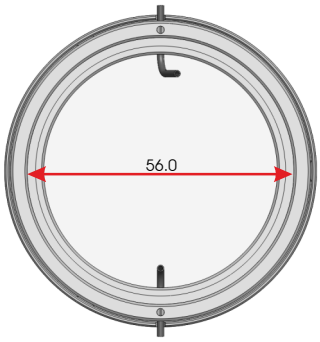dimensions in mm

- For "60" Petri dishes with an outer diameter:
 Dmax = 56.0 mm
 Dmin = 50.5 mm

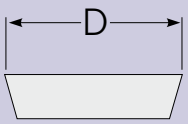
"60" Petri dish without cover


Assembly

- Set the sterilized perfusion ring onto the Petri dish "60"
- Set the cover with glas insert onto the perfusion ring
- Plug the Tygon® tubes onto the canal tubes

Application example

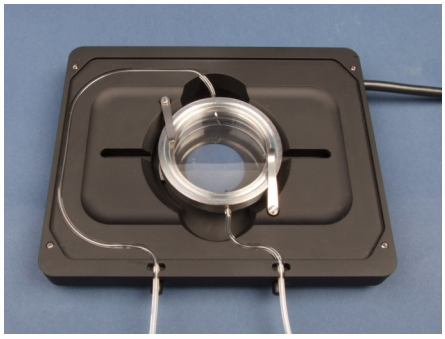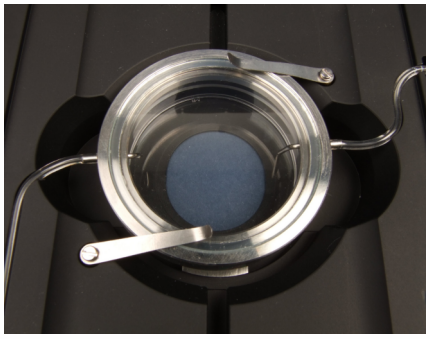- Petri dish "60" with perfusion ring and cover insered into the Heating Insert P 2000

This Perfusion-Set is compatible to a large variety "60" Petri dishes.
Please do not hesitate to contact us regarding a cover with perfusion adapter for other Petri dishes!

Version 1.2, 03-2019
www.pecon.biz

PECON®
Systems for cell culture and microscopy

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